

KIT

ALGLASS

ALUMINA GLASS COMPOSITE

FOR STRUCTURE AND COPING MAKING
NOW AVAILABLE IN SHORT KIT VERSION



M.S. 1180064100001



METAL FREE IN YOUR HAND



Alumina sintered at 1700°C is one of the most resistant materials known nowadays. However, the high temperature and thermal contraction inherent to the process make the use of this pure material impracticable for dental structure making.

Materials usually called glasses have low fusing point but, on the other hand, have small mechanical resistance. The glass infiltrated in a pre-sintered alumina matrix results in a composite that fulfill the mechanical requirements and biocompatibility in a "Metal-Free" technique for structures and coping making.

Material and Chemical Engineers Researchers from Brazilian Universities together with specialists in Dental Prosthesis and Celmat ® developed a family of products that uses this technology, making this idea available to the market at low cost, high quality and reliability.



PLASTER

Kit - 2 package - 200g



ALUMINA

Kit - 10 packages of Alglass Alumina containing 38g
Short Kit - 03 packages of Alglass Alumina containing 38g



LIQUID

Kit - 10 bottles of Alglass liquid containing 5 ml each
Short Kit - 03 bottles of Alglass liquid containing 5 ml each



DISPERSANT

Kit - 1 bottle of Alglass dispersant containing 10 ml
Short Kit - 1 bottle of Alglass dispersant containing 10 ml




CRACK CHECKER

Kit - 1 bottle of Alglass crack checker containing 10 ml
Short Kit - 1 bottle of Alglass crack checker containing 10 ml



BEAKER

Kit - 1 Beaker of 50 ml
Short Kit - 1 Beaker of 50 ml

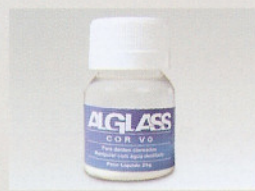

GLASS "V0" COLOR

Kit - 1 bottle of Alglass glass "V0" color containing 25g
Short Kit - bottle of Alglass glass "V0" color containing 17g



GLASS "V1" COLOR

Kit - 1 bottle of Alglass glass "V1" color containing 25g
Short Kit - bottle of Alglass glass "V1" color containing 17g



GLASS "V2" COLOR

Kit - 1 bottle of Alglass glass "V2" color containing 25g
Short Kit - bottle of Alglass glass "V2" color containing 17g



BRUSH

Kit - 1 brush
Short Kit - 1 brush



SPATULA

Kit - 1 Glass Spatula
Short Kit - 1 Glass Spatula



SYRINGE

Kit - 1 Syringe of 20 ml
Short Kit - 1 Syringe of 20 ml

CELMAT